


AirMini P10 Starter Pack

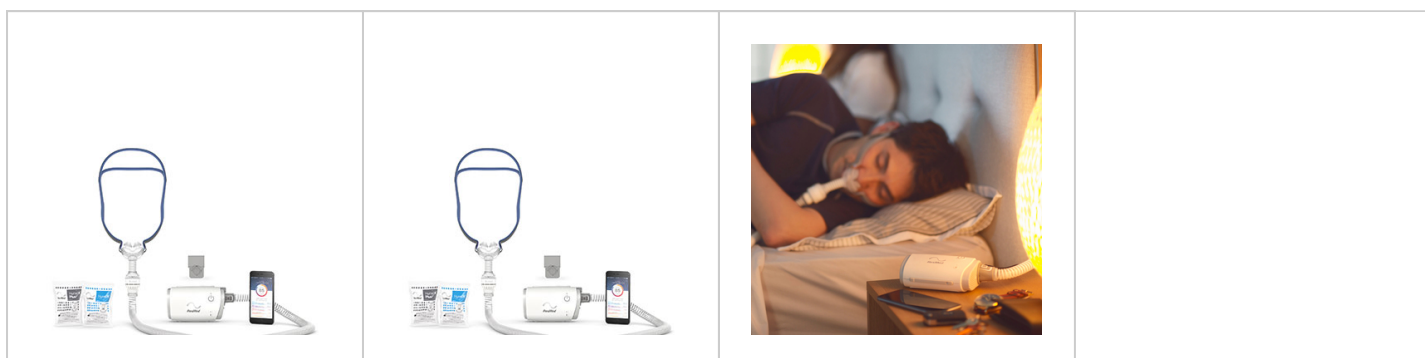

The P10 starter kit includes the AirMini Machine, P10 Set up pack (tubing connection), Mount system and P10 mask and includes the small, medium and large pillows. Small in size and smart with features the ResMed AirMini provides you with the quality sleep apnea therapy you've come to expect from ResMed- only in a single, complete package that's compact and lifestyle-friendly. Whether you are looking for a portable travel CPAP or a primary device for your sleep apnea treatment, the AirMini could be the solution you've been looking for!

The best thing of all is you can buy and then try and if you are unhappy with the device we have a **30 Day money Back Guarantee** - no questions asked.

This is the best possible solution for frequent travellers and business men who don't want to have to pack their bulky machine all the time. It provides all the comforts of your standard machine and can be used as a front line CPAP every night

When ordering this machine it comes as a Bedside starter kit which includes the AirMini device, mask starter pack, mask specified and bedside mounting kit. It is very important to order the correct mask as the kit is specific to the mask. The AirMini will only operate with the latest Range of Resmed masks and that is why buying as a kit you need to specify the mask. You will not be able to use your current mask if you have one already as the connection to the tube will not be compatible.

Key Features

It's small! The AirMini is the smallest CPAP machine in the world*. It measures 13.6 cm wide and weighs only 300 grams!

It's innovative! ResMed have developed a tube and connector with a built-in vent. To deliver quality performance at such a small size, these connectors will only work with ResMed's AirFit™ 20 series of full face masks or nasal (N20) masks, as well as with the ResMed AirFit P10 Pillow mask

It's feature packed! This compact AirMini CPAP machine includes CPAP (fixed pressure), APAP (automatic) and APAP for Her modes all in the same device. It also offers features such as Auto Ramp,

Expiratory Pressure Relief plus Smart Start and Stop. And it can provide a maximum pressure of up to 20!

It's waterless! The AirMini uses heat moisture exchange (called HumidX) to deliver waterless humidification.

It's smart! You can conveniently control the AirMini via a smartphone or tablet app which also allows you to monitor your therapy.

It has a 30 day money back guarantee! AirMini delivers the same quality of sleep therapy that you've grown to trust from ResMed. That's why ResMed are offering you a 30 day money back guarantee if you are not completely satisfied..

Product Information

Frequently Asked Questions

What is the AirMini's maximum pressure?

The AirMini provides a pressure range from 4-20cm H₂O

How do I set up the AirMini with my therapy settings?

We recommend setting up a time to call MyCPAP. We will then set the machine up using your therapy settings and discuss how it works. If you are purchasing via our online store, you can select your preferred collection location when you go to the check-out page.

What modes does the ResMed AirMini machine come with?

The AirMini machine comes with all three CPAP modes: AutoSet, AutoSet for Her and CPAP. It uses the same algorithms used in ResMed's AirSense 10 machines.

Why is there no screen?

The AirMini has been designed to support both simplicity and a richer experience depending on your preference. To support this, ResMed designed AirMini with a simple one-button interaction for users who prefer simplicity and an accompanying smartphone app to unlock a richer, more engaging therapy experience.

Do I need a smartphone to use the AirMini?

No. A smartphone is only needed to set up the AirMini, not to use it. Once your AirMini has been set up with your therapy settings at one of our stores, there is no need for a smartphone to use the AirMini. You can just use the single on/off button to operate the AirMini.

If I use the AirMini in other countries, will the power supply unit automatically convert?

Yes, the design of ResMed's AirMini machine accepts an AC input range of 100-240V and a DC output of 24V. However, purchasing the correct AirMini power supply unit for the country you are visiting is recommended to ensure correct performance of your AirMini machine.

Is there a battery available for the AirMini?

At this stage a battery isn't available for the AirMini. However, ResMed are working to have an externally battery available later in 2017.

Is there a DC-DC converter available?

This isn't currently available, but will be available soon.

What is ActiveAir technology?

ActiveAir technology is a mask vent that uses a silicone membrane to regulate vent flow at all treatment pressures.

How loud is the ResMed AirMini?

The sound pressure level for the AirMini is 29dB

Can I use the AirMini machine when I travel on an airplane?

If you have access to power, you can use your AirMini machine when you travel on an airplane. Please set it to airplane mode before you travel (refer to travelling section of your AirMini user guide for instructions on entering and exiting airplane mode).

How much does the AirMini cost?

The RRP, set by ResMed, is \$2100 AUD. Check out MyCPAP ! This is for the starter kit which includes all items, such as Mask, Tubing, Power adapter, that you will need to start using the AirMini straight away

Why can I not use my current mask with AirMini?

The ActiveAir and HumidX technology is a new technology only implemented in the AirMini. As a result the AirFit 20 series masks and AirFit P10 for AirMini are the only masks designed to work with the new AirMini system.

Is the AirMini compatible with nasal pillows made by other companies?

The AirMini will only work with ResMed AirFit F20, N20 or P10 CPAP masks

Does the AirMini have humidification?

AirMini uses atmospheric moisture as well as moisture from your breath to humidify the airflow. This waterless humidification technology is known as heat moisture exchange (HumidX).

What is HumidX?

HumidX is the name ResMed has given the AirMini humidification. It utilises a disposable capsule that sits in the tube system which is located right before the air enters the mask. Capsules are available in two strengths: HumidX and HumidX Plus.

Which masks can I use HumidX with?

You can use HumidX with ResMed's Nasal (AirFit N20 and N20 for Her) and Nasal Pillows (AirFit P10). HumidX is not currently available for Full Face masks.

How often does the HumidX need to be replaced?

Replace the HumidX capsules 30 days after the packaging is opened. If the HumidX was not used every

day during the 30-day period, it must still be replaced.

Can I clean the HumidX and re-use it?

No. The HumidX cannot be cleaned or washed and must be replaced 30 days after opening the packaging.